

# Creation of output spreadsheet

Script functions on the premise that input files have 5 sheets organized as: Result, Controls, Calibrants, Rbl, Method

Start with sample results

Extract the last run of each sample (the run with greatest dilution ratio)

Check how many duplicates were run for each sample (same SampleID, position, dilution)

Script makes sure that each run has 2 duplicates. If 2<sup>nd</sup> duplicate is missing, script creates a dummy duplicate (concentration = NA)  
If more than 2 duplicates were run, script takes note and later flags as a “Multi-run sample”

Assign additional flags to each sample run (see Final Flags)

Use the concentrations from the 2 duplicates to calculate

- mean
- difference and percent difference
- standard deviation and relative standard deviation


If either of the concentrations < 0, the mean is set to -111

Assign flags to the newly calculated values (see Final Flags)

If only 1 duplicate was run, the calculated difference, percent difference, standard deviation, and relative standard deviation are NA

Place all data together in one final output spreadsheet

# Creation of output plots

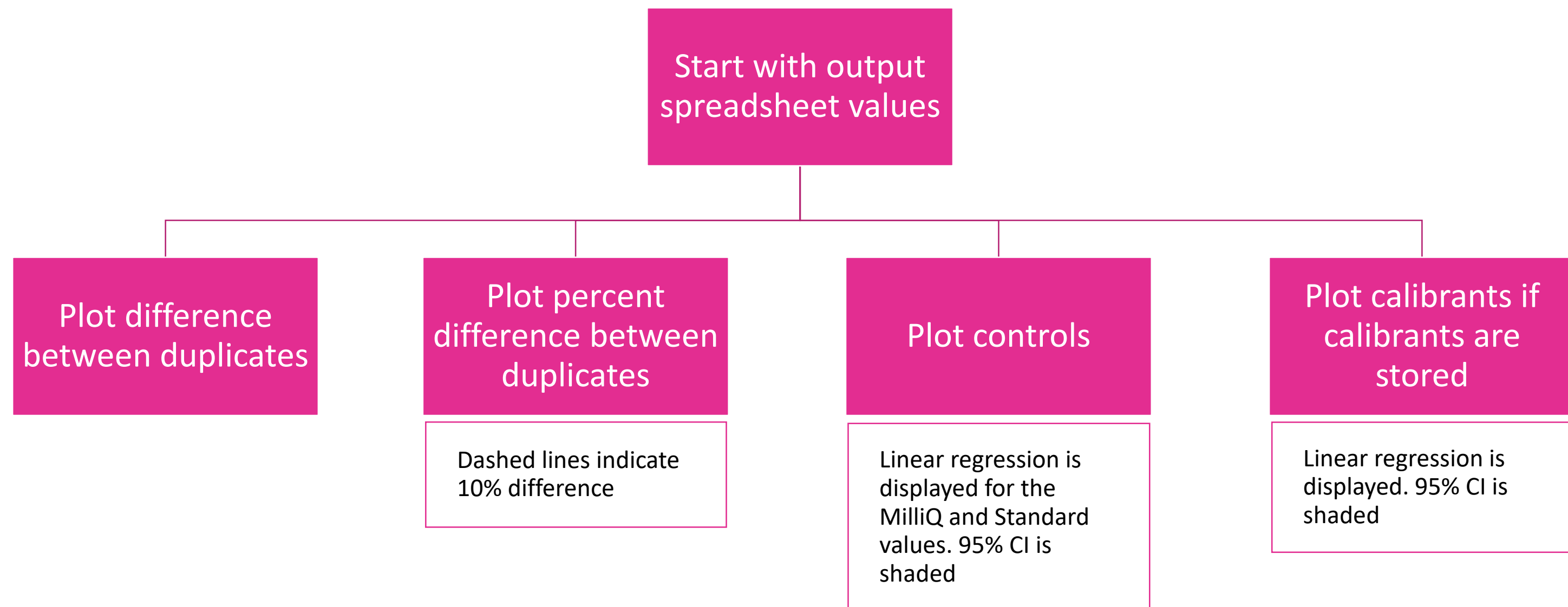

# Conversion of SC flags

SC flag	Expanded flag
LL	Linearity Low
LH	Linearity High
EPL	End-Point Limit
INV	Inversion
><	Out of Calibration Curve
SS	Short on Sample
SR	Short on Reagent
H	High

# Final Flags

## Flags for individual runs (columns: Flags Dup x)


Flag	Meaning
Linearity Low	Abs value is low and outside of the calibration curve
Linearity High	Abs value is high and outside of the calibration curve
End-Point Limit	Sample sat longer than recommended while being read
Inversion	Negative absorbance
Out of Calibration Curve	Result plotted outside of the calibration curve
Short on Sample	Run could not complete due to insufficient sample
Short on Reagent	Run could not complete due to insufficient reagent
High	Sample was above the calibration curve
Run was not completed	Will apply only to _Dup rows. Indicates that _Dup row is missing in input
< 0 (actual value: xxx)	A negative concentration was recorded. Concentration value is stored in flag
< MDL	A concentration < MDL was recorded

## Flags for comparison of duplicates (column: Flags)

Flag	Meaning
Perc difference > 10 & both conc above MDL	Percent difference is violated and both concentrations >= MDL
Perc difference > 10 & one or both conc below MDL	Percent difference is violated and one or both concentrations < MDL
Both conc below MDL	Percent difference not violated and concentrations < MDL
Rerun dilution	SC automatically reran this sample with a different dilution
Multi-run sample	User reran this sample with the same dilution

# Callstack

For those who need to go into the code

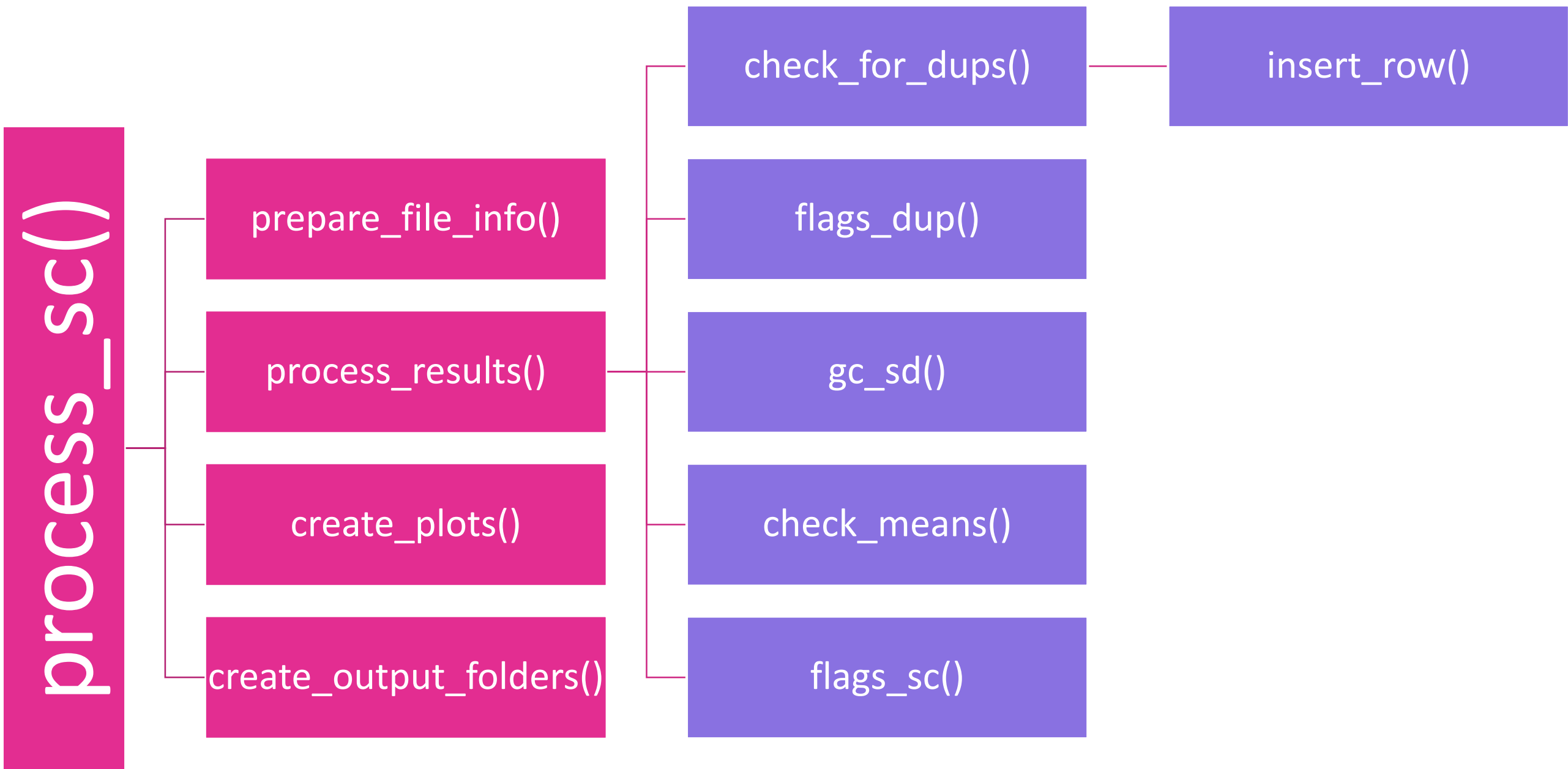

function defined in  
sc\_functions.R

function defined in  
sc\_helper\_functions.R